

Urimeter



The sterile, accurate
closed collection system

flexicare
TOTAL QUALITY - TOTAL CARE

Urimeter

The Flexicare Urimeter is a sterile closed collection system with an integral measuring chamber. The vertical drainage system and incorporated non return valves (NRV's) ensure the measuring chamber empties fully once the reading has been taken to ensure accurate measurements especially with low urine output.




Anti-kink soft tubing
Anti-kink attachment allows unimpaired flow whilst our soft tubing aids patient comfort.

Hydrophobic filters
Allowing for any air build up to escape without allowing any proteins or fluids to pass through.



Need-free sample port
Eliminates the risk of needle stick injury.

Vented valve system
Maintains urine flow whilst reducing hydrostatic pressure, minimising damage to the bladder wall and reducing catheter blockage.

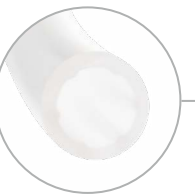

Hexagonal tubing
Reducing surface tension allowing urine flow and increased accuracy.



Secure fastening
The integral hangers and detachable straps ensure the Urimeter is securely attached to any bed frame.



Easy to read chamber
The three chamber system allows accurate readings of small volumes of urine with a smooth action 180° tap.

Safe and Secure Drainage
Non-return valves reduce the risk of ascending infection and allows low pressure drainage without compromising the anti-flux property.



90° Lever Tap
90° lever tap for easy drainage with protective housing to reduce contamination and accidental opening.



Vertical Drainage
The vertical drainage ensures complete emptying of the measuring chamber for accurate urine readings.

Clear print and date box
Clear indications of use and a date box allow information to be recorded and length of use monitored.

DESCRIPTION	CODE	QTY
Urimeter with 150mm Tubing	00-1227	20