

# Urimeter



The sterile, accurate  
closed collection system

**flexicare**  
TOTAL QUALITY - TOTAL CARE

# Urimeter

The Flexicare Urimeter is a sterile closed collection system with an integral measuring chamber. The vertical drainage system and incorporated non return valves (NRV's) ensure the measuring chamber empties fully once the reading has been taken to ensure accurate measurements especially with low urine output.




**Anti-kink soft tubing**  
Anti-kink attachment allows unimpaired flow whilst our soft tubing aids patient comfort.

**Hydrophobic filters**  
Allowing for any air build up to escape without allowing any proteins or fluids to pass through.



**Need-free sample port**  
Eliminates the risk of needle stick injury.

**Vented valve system**  
Maintains urine flow whilst reducing hydrostatic pressure, minimising damage to the bladder wall and reducing catheter blockage.

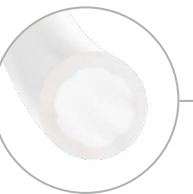

**Hexagonal tubing**  
Reducing surface tension allowing urine flow and increased accuracy.



**Secure fastening**  
The integral hangers and detachable straps ensure the Urimeter is securely attached to any bed frame.



**Easy to read chamber**  
The three chamber system allows accurate readings of small volumes of urine with a smooth action 180° tap.

**Safe and Secure Drainage**  
Non-return valves reduce the risk of ascending infection and allows low pressure drainage without compromising the anti-flux property.



**90° Lever Tap**  
90° lever tap for easy drainage with protective housing to reduce contamination and accidental opening.



**Vertical Drainage**  
The vertical drainage ensures complete emptying of the measuring chamber for accurate urine readings.

**Clear print and date box**  
Clear indications of use and a date box allow information to be recorded and length of use monitored.

DESCRIPTION	CODE	QTY
Urimeter with 150mm Tubing	00-1227	20